

# MSCs-derived exosomes

REF	SIZE	DESCRIPTION
ExoMSC	request information	Adipose-derived Mesenchymal stem/stromal cells (MSCs) derived exosomes.

## 1. INTRODUCTION

Exosomes are small extracellular vesicles that are released from cells upon fusion of an intermediate endocytic compartment, the multivesicular body (MVB) <sup>(1)</sup>, with the plasma membrane. They are thought to provide a means of intercellular communication <sup>(2,3)</sup> and of transmission of macromolecules between cells allowing the spread of proteins, lipids, mRNA, miRNA and DNA and as contributing factors in the development of several diseases.

Exosomes can also modulate cancer microenvironment <sup>(4)</sup> and the immune response <sup>(5,6)</sup>. In this sense, it is well known that Mesenchymal stroma/stem cells (MSCs) are clinically useful for cell-based therapy, however some studies have shown that MSC-derived Exosomes emulate the effect of MSCs in various experimental models, stimulating cell proliferation and repair <sup>(7)</sup>.

## 2. PRODUCT DESCRIPTION

Lyophilized/Frozen exosomes (~1x10<sup>12</sup>) from adipose-derived Mesenchymal stem/stromal cells (MSCs) <sup>(8,9)</sup>. Exosomes are isolated by differential ultracentrifugation <sup>(10)</sup>.

- **Tested application:** In-vitro functional assay, Flow Cytometry (FMC), Nanoparticles Tracking Analysis (NTA, Nanosight), Western Blot (WB), BCA Protein Assay.
- **Species reactivity:** Human
- **Presentation:** Lyophilized/Frozen
- **Reconstitution of lyophilized Exosomes:** For reconstitution, we recommended adding sterile, distilled water to achieve a final exosome concentration of 1µg/µl (e.g., for 100 µg standard, add 100 µl of dH<sub>2</sub>O). After the addition of water, recap vial and briefly vortex making sure that the liquid has been gently distributed and has covered the entire inside of the vial.  
After vortexing, make sure that the solution is collected at the bottom of the vial, if not, centrifuge shortly the vial solution. Now the standard is ready to use.

## 3. APPROPRIATE STORAGE AND HANDLING CONDITIONS

Lyophilized exosomes can be stored between 2°C and 8°C for up to 2 years without functional compromise. Frozen exosomes can be stored at -20 °C.

Immunostep recommends storing small, single –use aliquots of reconstituted exosomes, at – 20°C for up to one month or at – 80°C for longer periods, preferably in locations in frost-free freezers, without appreciable temperature fluctuation. This will minimize protein denaturation that can occur after multiple freeze/thaw cycles.

Reconstituted exosomes, store properly, are functionally guaranteed for up to six months from date of reconstitution.

Any unfrozen and/or unused exosome standard can be stored at 4°C for short term use (<1 week), and should not be re-frozen.

## 4. EVIDENCE OF DETERIORATION

Reagents should not be used if any evidence of deterioration is observed. For more information, please contact our technical service: [tech@immunostep.com](mailto:tech@immunostep.com)

## 5. WARRANTY

Warranted only to conform to the quantity and contents stated on the label or in the product labelling at the time of delivery to the customer. Immunostep disclaims hereby other warranties.

Immunostep's sole liability is limited to either the replacement of the products or refund of the purchase price.

## 6. PERFORMANCE DATA

At the time of exosome collection, for each batch, MSCs phenotype analysis is performed by FCM in order, to guarantee the right signature described by The Mesenchymal and Tissue Stem Cell Committee of the International Society for Cellular Therapy (ISCT) <sup>(11)</sup>. Thus MSC must express CD105, CD73 and CD90, and lack expression of CD45, CD34, CD14 or CD11b, CD79alpha or CD19 and HLA-DR surface molecules (Fig. 1).


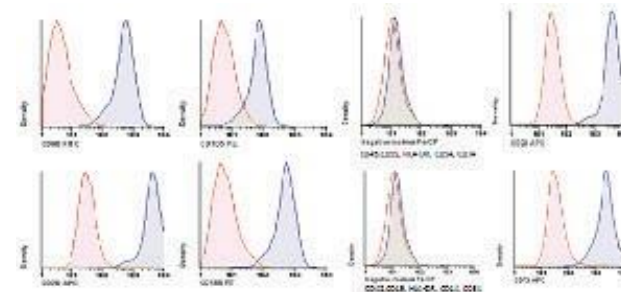


Figure 1: The following histograms correspond to the analysis of antigen expression in an expansion of Adipose-derived Mesenchymal stem/stromal cells MSCs (passage 2) cultured for exosome production.

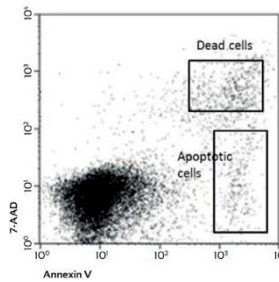

The percentage of apoptotic and necrotic MSCs at the time of collection is also analyzed by FCM (Fig.2) to verify that exosome preparation is not cross-contaminated with apoptotic bodies.

Exosome batches are checked and compared for the presence of the CD63 and CD9, a common exosome marker, by FCM and WB (Fig.3).

Figure 2: Apoptotic and necrotic MSCs percentage analysis by FCM at the time of supernatant collection.

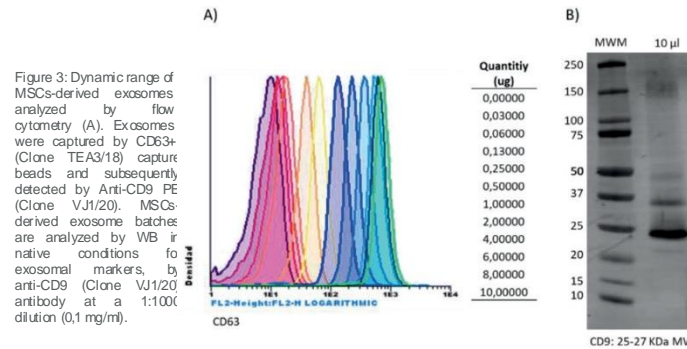


Figure 3: Dynamic range of MSCs-derived exosomes analyzed by flow cytometry (A). Exosomes were captured by CD63+ (Clone TEA3/18) capture beads and subsequently detected by Anti-CD9 PE (Clone VJ1/20). MSCs-derived exosome batches are analyzed by WB in native conditions for exosomal markers, by anti-CD9 (Clone VJ1/20) antibody at a 1:1000 dilution (0,1 mg/ml).

All exosome batches are also subjected to NTA Analysis to test for particle size and concentration (Fig.4).

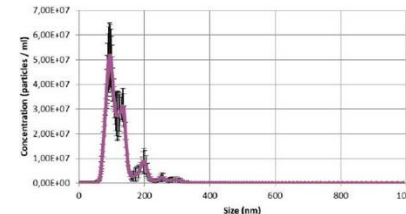


Figure 4: MSCs-derived exosome analysis for particle size and concentration by NTA, NanoSight LM10HSB. Analysis is carried out with 1 µl of purified exosomes diluted in 999 µl of HEPES buffer (dilution 1:1000). The purified exosomes showed a size distribution profiles, with peak diameters from 50 – 150 nm and concentrations about 1x10<sup>12</sup> exosomes/ml.

Finally, exosomes derived from MSCs are functionally analyzed in an in vitro assay (Fig. 5), in which uptake cell and internalization exosomes are checked on wide range of different cell cultures.

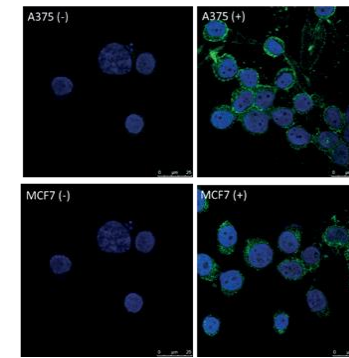


Figure 5: Cellular uptake of exosomes. MSCs-derived exosomes are isolated and labeled with PKH67 and DAPI. The labeled exosomes are cultured with several cell lines (here are shown A375 and MCF7) for 17 hours and analyzed by confocal microscopy. PKH67 bright spots were seen in the cells cultured with loaded exosomes, and not seen in control cells, whilst DAPI stains nucleic acid in both.

## 7. ADDITIONAL INFORMATION

For research use only. Not for diagnostic use.

Not for resale. Immunostep will not be responsible of violations that may occur with the use of this product. Any use of this product other than the specified in this document is strictly prohibited.

Unless otherwise indicated by Immunostep by written authorization, this product is intended for research only and is not to be used for any other purpose, including without limitation, for human or animal diagnostic, therapeutic or commercial purposes.

Please, refer to [www.immunostep.com](http://www.immunostep.com) technical support for more information.

## 8. EXPLANATION OF SYMBOLS

	Fluorochrome
	Product reference
	Content for <n> analysis
	Regulatory Status
	Research Use Only
	Manufacturer

## 9. REFERENCES

1. Yáñez-Mó M, Siljander P, Andreu Z, Bedina Zavec A, Borràs F, Buzas E et al. Biological properties of extracellular vesicles and their physiological functions. *Journal of Extracellular Vesicles*. 2015;4 (1):27066.
2. Pitt JM, André F, Amigorena S, Soria JC, Eggermont A, Kroemer G, Zitvogel L. Dendritic cell-derived exosomes for cancer therapy. *J Clin Invest*. 2016.
3. Tkach M, Théry C. Communication by Extracellular Vesicles: Where We Are and Where We Need to Go. 2016 *Cell* 10;164(6):1226-32.
4. Becker A, Thakur BK, Weiss JM1, Kim HS, Peinado H, Lyden D Extracellular Vesicles in Cancer: Cell-to-Cell Mediators of Metastasis. *Cancer Cell* 2016 Dec 12;30(6):836-848.
5. López-Cobo S, Campos-Silva C, Valés-Gómez M. Glycosyl-Phosphatidyl-Inositol (GPI)-Anchors and Metalloproteases: Their Roles in the Regulation of Exosome Composition and NKG2D-Mediated Immune Recognition. *Front Cell Dev Biol*. 2016 Sep 12;4:97.
6. Jonathan M. Pitt, Guido Kroemer, Laurence Zitvogel Extracellular vesicles: masters of intercellular communication and potential clinical interventions. 2016 *J Clin Invest* 2016;126(4):1139-1143
7. Vishnubhatla I, Corteling R, Stevanato L, Hicks C, Sinden J The Development of Stem Cell-Derived Exosomes as a Cell-Free Regenerative Medicine. *Journal of Circulating Biomarkers*. 2014;3:2.
8. Campos S, Suárez H, Jara-Acevedo R, Linares-Espinós E, Martínez-Piñero L, Yáñez-Mó M, Valés-Gómez M. High sensitivity detection of extracellular vesicles immune-captured from urine by conventional flow cytometry. *Sci Rep*. 2019; Feb 14;9(1):2042.
9. Jara-Acevedo R, Campos-Silva C, Valés-Gómez M, Yáñez-Mó M, Suárez H, Fuentes M. Exosome beads array for multiplexed phenotyping in cancer. *J Proteomics*. 2019; Apr 30;198:87-97.
10. Théry C, Amigorena S, Raposo G, Clayton A. Isolation and characterization of exosomes from cell culture supernatants and biological fluids. *Curr Protoc Cell Biol*. 2006 Apr;Chapter 3:Unit 3.22.
11. Dominici M, Le Blanc K, Mueller I, SlaperCortenbach I, Marini F, Krause D, et al. Minimal criteria for defining multipotent mesenchymal stromal cells. The International Society for Cellular Therapy position statement. *Cytotherapy*2006;8(4):315-7

## 10. MANUFACTURED BY



### IMMUNOSTEP S.L.

Address: Avda. Universidad de Coimbra, s/n  
 Cáncer Research Center (C.I.C)  
 Campus de Unamuno  
 37007 Salamanca (Spain)  
 Telf./fax: (+34) 923 294 827  
 E-mail: [info@immunostep.com](mailto:info@immunostep.com)  
[www.immunostep.com](http://www.immunostep.com)